

No Cleaning — No Cross Contamination — Reduce Costs

Jumbo Pipettes


Low Cost Disposable Liquid Sampler

Jumbo Pipettes are ideal for sampling low viscosity liquids (such as water). They have tapered ends to allow a large volume sample to be taken. The ends have small holes in to ensure that the operator has complete control over the product flow when discharging the sampler. Also minimises/eliminates dripping.


- Simple to use
- Two sizes available
- Individually wrapped
- Assembled in a medical cleanroom
- Available pre-sterilised
- Made from HDPE (FDA acceptable grade - conforms to CFR 177.1520). The blue colouring is an FDA acceptable dye.
- Fully welded construction - no adhesives used
- ATEX rated for zones 0,1 & 2



Hold the Blue end - place thumb over hole to prevent the sample from running out

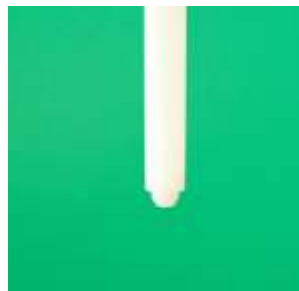

The white end goes into the product

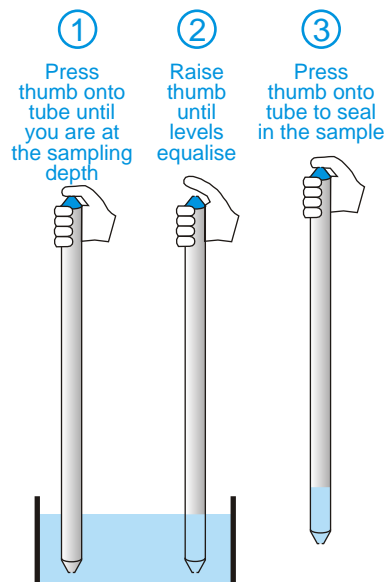

Jumbo Pipette
Straight sides
Tapered Ends

Jumbo Pipettes - large volume, with extra control



Part No.	Overall Length	Overall Diameter	Internal Diameter	Pre-Sterilised	Material	Max. sample Volume*	Quantity per Box
8051H-450	450mm	21mm	17mm	No	HDPE	80ml*	50
8051H-450S	450mm	21mm	17mm	Yes	HDPE	80ml*	50
8051H-900	900mm	21mm	17mm	No	HDPE	180ml*	50
8051H-900S	900mm	21mm	17mm	Yes	HDPE	180ml*	50

* The maximum sample volume will depend upon how deeply the sampler is inserted into the product